

## Puritan Bennett Outlet Upgrade Kits For PB Quick-Connect and DISS

### SPECIFICATION

**Outlet Upgrade Kit**

The Outlet Upgrade Kits shall permit the customer to update older series Puritan Bennett outlets to the new Series B configuration. The Outlet Upgrade Kits shall permit the customer to update early DISS or quick-connect outlets to either current DISS or quick-connect. Configurations that can be updated include the Litharged series manufactured from 1953 to 1963, the series 73 manufactured from 1963 to 1973, and the series A outlets manufactured from 1973 to 1988.

**Outlet Upgrade Kit Design**


Each upgrade kit shall be gas specific and require a connector kit that is specific to the outlet being replaced. In addition to the upgrade kit and proper connector, a standoff and cover plate shall be required to convert each outlet.

Item	Standard Console		
	Gas Service	DISS Upgrade Kit <sup>1,2,3</sup>	Puritan Bennett Quick-Connect Upgrade Kit <sup>1,2,3</sup>
U.S. Labeling	O <sub>2</sub>	<input type="checkbox"/> 6-232200-00	<input type="checkbox"/> 6-232020-00
	N <sub>2</sub> O	<input type="checkbox"/> 6-232201-00	<input type="checkbox"/> 6-232021-00
	AIR	<input type="checkbox"/> 6-232202-00	<input type="checkbox"/> 6-232022-00
	VAC	<input type="checkbox"/> 6-232203-00	<input type="checkbox"/> 6-232023-00
	WAGD	<input type="checkbox"/> 6-232209-00	<input type="checkbox"/> 6-232024-00
International Labeling	O <sub>2</sub>	<input type="checkbox"/> 6-232210-00	<input type="checkbox"/> 6-232025-00
	AIR	<input type="checkbox"/> 6-232212-00	<input type="checkbox"/> 6-232027-00
	VAC	<input type="checkbox"/> 6-232213-00	<input type="checkbox"/> 6-232028-00
Miscellaneous	Slide Bracket Upgrade Kit	<input type="checkbox"/> 6-232029-00	
	Series A Connector Kit	<input type="checkbox"/> 6-232032-00	
	Series 73 Connector Kit	<input type="checkbox"/> 6-232033-00	
	Litharged Connector Kit	<input type="checkbox"/> 6-232034-00	
	Buried Outlet Adapter Kit	<input type="checkbox"/> 6-232031-00	
	5" Standoff (Series A)	<input type="checkbox"/> 6-315014-00	
	5" Cover Plate (Series A)	<input type="checkbox"/> 6-325149-00	
	3 ¼" Standoff (Series 73 and Litharged)	<input type="checkbox"/> 6-315013-00	
	3 ¼" Cover Plate (Series 73 and Litharged)	<input type="checkbox"/> 6-325160-00	
	Staked Connector Tool Kit <sup>4</sup> (Required for Series 73)	<input type="checkbox"/> 6-138024-00	

**Note:**

1. Before ordering, confirm that the outlet series to be upgraded matches the illustrations provided on page 2.
2. Series A conversions require 5" standoffs and cover plates.
3. Series 73 and Litharged models require 3 ¼" standoffs and cover plates.
4. The staked Connector Tool Kit (required for Series 73) includes:  
 A 5/8"-24 bottom thread tap, a connector socket, and a 5/16" open-end wrench.

**Old Series Puritan Bennett Outlets**

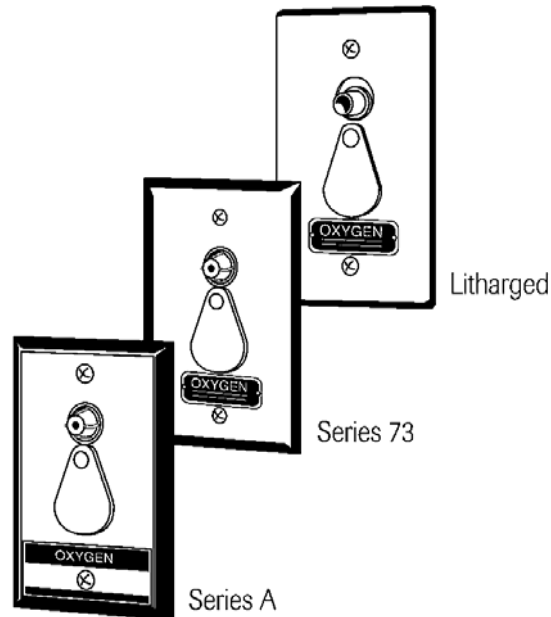

**Buried Outlet Adapter Kit (Optional)**

**NOTE:**

Some outlets may have been installed below the finished wall surface too far to permit upgrade installation. A modified adapter nut assembly, designed to compensate for this variance, is available (P/N 6-232031-00).

To determine if the modified adapter nut assembly is required, measure the distance from the end of the upgrade connector to the finished wall surface. If the dimension exceeds 13/16", the modified adapter nut assembly is required.

